



PITT-EASY®
Product Catalog



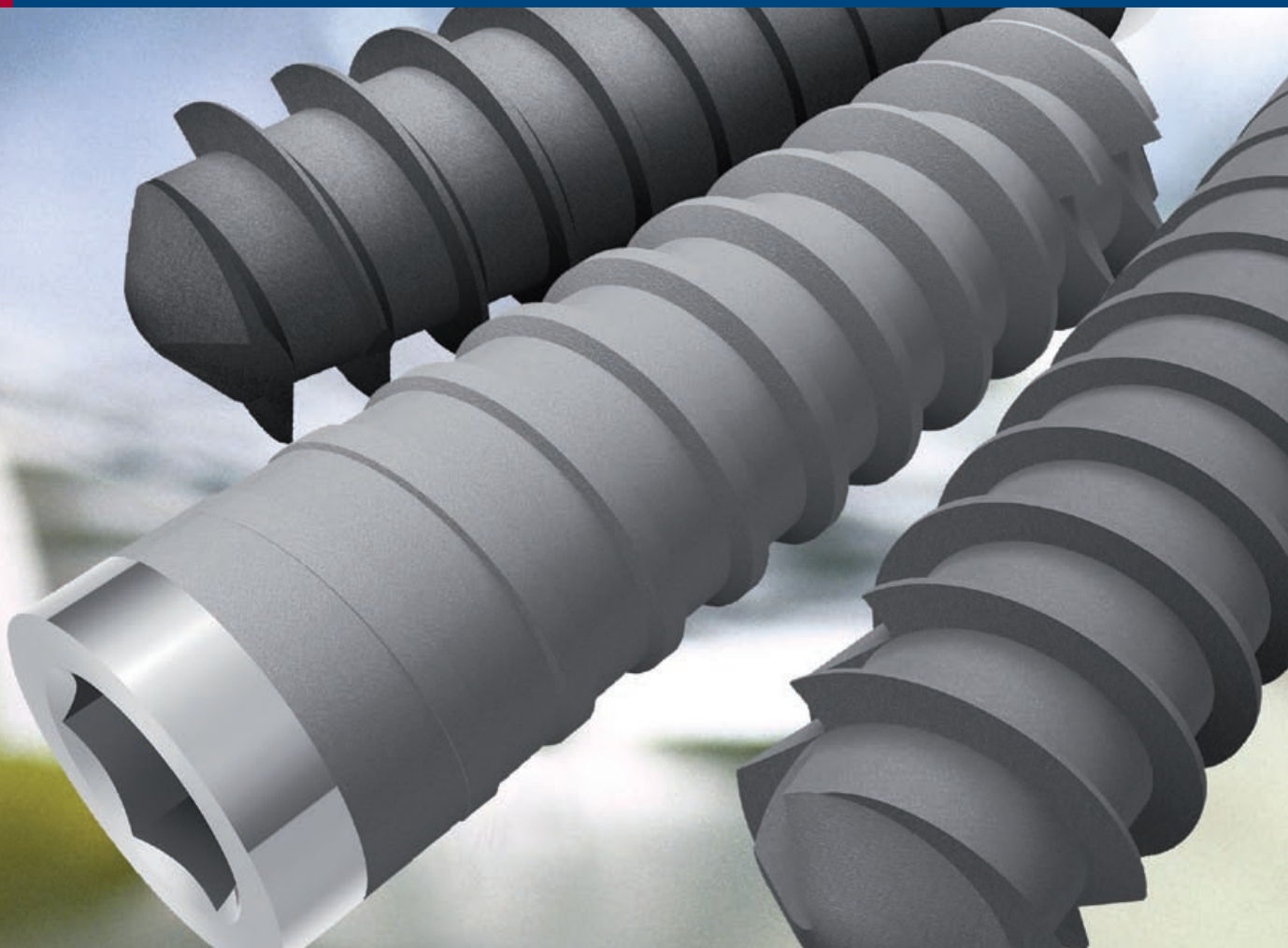

Contents

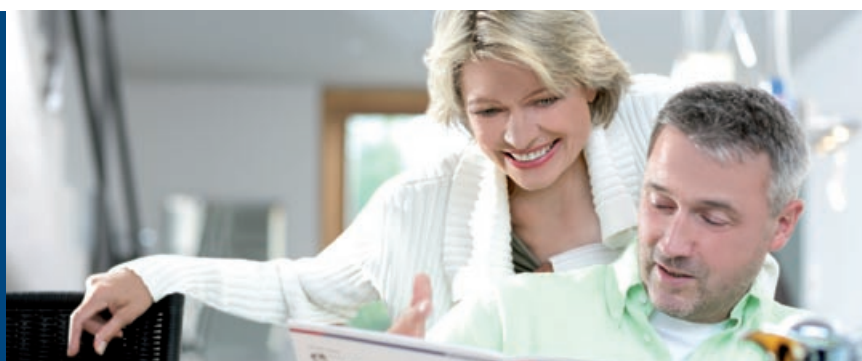

Implant system

- 4 Quality certified
- 5 Three surfaces at your choice
- 6 Color Guide System

Surgery

- 7 Pre-surgery planning
- 8 Implants with V-TPS, FBR and PURETEX®-Surface
- 9 Simultaneous Drills
- 10 Initial Drills
- 11 Universal Drills
- 13 Surgical Trays
- 16 Instruments and Accessories

Prosthetics

- 18 Impression
- 20 Abutment Trial Set
- 21 A.G.T. Line
- 22 V.D.L. Esthetic Line
- 24 Castable Abutment System
- 26 CAD/CAM-Base
- 28 Individual Abutments
- 30 Paracentric Line
- 32 Ball Head Abutment System
- 34 Locator Abutment System
- 35 Standard Abutment with hex
- 36 Prosthetic Accessories
- 38 Table Prosthetic Lines and Abutments
- 40 Prosthetic Kit
- 41 Table Prosthetic Drivers and Screws
- 42 Prosthetic Drivers, Instruments and Accessories

Safety Data Sheets

The material composition of all products is documented in the form of Safety Data Sheets. These are available upon request.

Internet

All print material released by Sybron Implant Solutions is available online, too. Product catalogs and insertion brochures as well as information for patients or current case reports – it takes you only 4 steps to get to the print product desired. Just click under www.sybronimplants.de

Implant system

Quality certified

Our quality assurance and control system has been certified under DIN EN ISO 13485:2003. This standard provides for the most exacting quality management requirements to be implemented and maintained at all levels concerned. At the same time, the conformity with the requirements of Directive 93/42/EEC concerning medical products was confirmed and all our products show the CE mark.

Each implant is individually sealed in a gamma-sterilized double safety package and sterility is assured for five years.

The package contains instructions for use with a guidance for insertion to ensure safe application, as well as detachable labels for the implant passport or the patient files.

The labels and the outside packaging contain all the important information: implant length, implant diameter, sterility expiration, item number and lot number (enables tracking of each implant back to its respective lot of the titanium manufacturer). In addition to the implant packages, also the boxes of

the drills, prosthetic abutments, drivers and model analogs contain two copies of the corresponding label.

The implant is kept in its acrylic vial until insertion. The color-coded screw cap contains a color-coded implant cover screw and indicates diameter and length of the implant.

For insertion, the implant is removed from the acrylic vial together with the screw cap, and is inserted into the prepared bone bed by two to three rotations. This procedure prevents contamination of the implant body.



Three surfaces at your choice

PURETEX® surface

PURETEX®, the biomimetic implant surface

Puretex has opened a completely new spectrum for implant therapy. Puretex derives from tex = texture which is the surface topography and from pure for purity. We have implemented a new decontamination technology which creates an implant surface of utmost purity.


The comprehensive studies and test series in close cooperation with university institutes have resulted in the

V-TPS (Vacuum-Titanium Plasma Spray) surface


development of our technology. We have succeeded in creating a biomimetic titanium surface by a purely subtractive process.

A clinical evaluation confirmed this excellent healing property. There was no horizontal bone resorption as with other implants.

FBR (Fast Bone Regeneration) surface

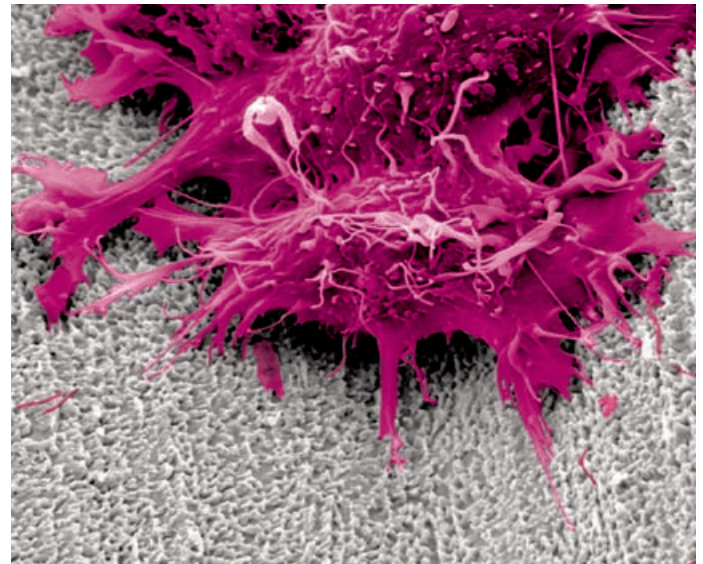

Advanced spreading, stage 4, formation of extremely numerous filopodiae with intensive contact to the biomimetic structure in immediate vicinity of the cell, complete flattening (magn. factor 1:5000, illustration is colorized)

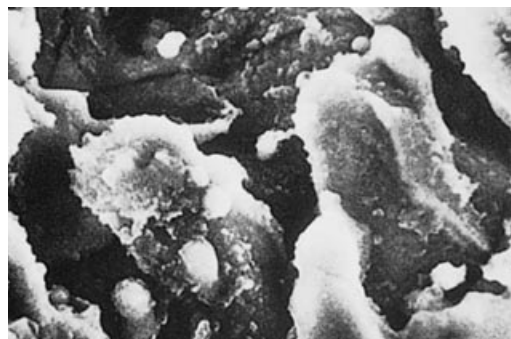

V-TPS surface (magn. factor 1:3000)

PITT-EASY® with V-TPS surface for approved therapy

The V-TPS surface enlarges the implant surface sixfold and ensures a contact osteogenesis and safe osseointegration. The special implant shape and the multiple diameters and lengths, developed for each indication, guarantee a quick, uncomplicated bone regeneration. The safe osseointegration has been documented for nearly 30 years. Advantage: The vacuum causes a homogeneous gap-free implant interface adherence with finest insoluble particle connection.


PITT-EASY® with bioactive FBR surface for early loading after six weeks

The FBR (Fast Bone Regeneration) surface consists of 100% calcium phosphate. It is additionally applied on the tried-and-tested fine granular V-TPS surface in a thin layer (15 – 20 µm). The FBR coating displays a marked capillary effect for an intensive, premature anchorage.

Pitt-Easy implants with FBR surface achieve a 23% higher rate of bone deposit than with other, rough-structured surfaces. This indicates early

loading is possible in a connection construction, i.e. 6 weeks post-op.

FBR – for patients who wish to shorten treatment schedules.

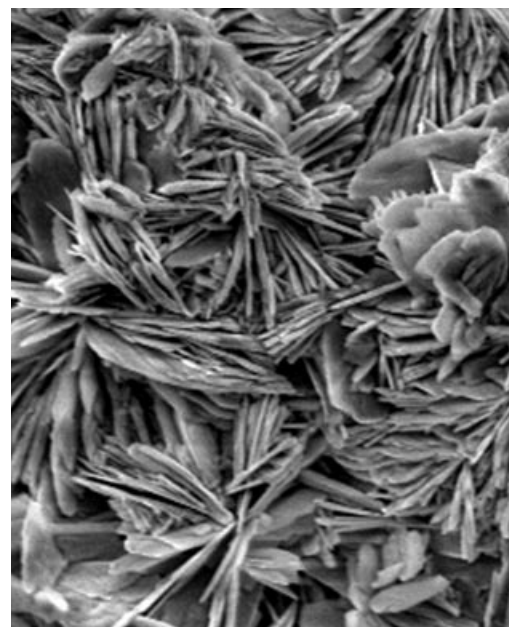


Fine crystalline surface structure FBR (magn. factor 1:2000)

Color Guide System

The Pitt-Easy implants are available in five diameters in order to meet the different alveolar ridge widths.

The abutments, impression copings and implant analogs are color-coded in diameters 3.75 mm (yellow) and 4.90 mm (blue). Healing abutments and Cover screws are color-coded for each diameter.

The surgical trays are also color-coded and facilitate the handling for dentists and assistants as well as for technicians.

Diameter	Color code	Implant color marking	Surgical tray marking
Ø 3.25 mm	green		
Ø 3.75 mm	yellow		
Ø 4.00 mm	red		
Ø 4.90 mm	blue		
Ø 6.50 mm	grey		

Legal notice

The implant may only be used by doctors who are suitably familiar with the system. The instructions (in this product catalog, the brochure "insertion technique", and the instructions for use) must be strictly complied

with. Continued training on this implant system is absolutely essential. We reserve the right to modify or improve the products in line with technical progress.



Surgery



Pre-surgical planning	Item No.	Ø mm
Drilling Guide Set set of 2 drilling guides each, Ø 3.25 – 6.50 mm, 5 positioning aid pins	90317	
Drilling Guides (package of 2) for the preparation of drilling templates (internal Ø 2.00 mm)	90306 90308 90310 90312 90314	Ø 3.25 Ø 3.75 Ø 4.00 Ø 4.90 Ø 6.50
Positioning Pins (package of 2) length 35 mm, for fixation of the drilling guides on the model	90316	Ø 2.00
Twist Drill Special for use with drilling guides, laser marked 8 – 20 mm, cutting at the apical portion, the remaining part of the cutting edge is dulled to avoid abrasion within the drilling guide	90009	Ø 2.00
X-Ray Indicator for selection of implant lengths and diameters, for placing on the panoramic radiograph, magnification factor 1:1/1:1.25/1:1.40	90074	X-Ray Indicator
X-Ray Balls package of 6	90150	Ø 5.00
Osteometer to determine the alveolar bone profile pre-surgical, for the preselection of the implant diameters	1210	



X-Ray Balls



| PITT-EASY® Implants

Implants

Item No.

Length mm

Pitt-Easy implants are delivered with a pertaining cover screw, stored in the lid of the acrylic vial. For reasons of hygiene, only single use of the Retaining screw is permissible.

The cover screws have an internal hex. The pertaining inbus driver is designed for a friction-fit.

When placing an order, please specify if you want implants with V-TPS surface (Item No. plus V-TPS), implants with FBR surface (Item No. plus FBR) and PURETEX surface (Item No. plus PUR). Please use the optimal implant length for each case.

Implants of diameter 3.25 mm are basically not suitable for single tooth replacement in the molar region.

Ø 3.25 mm PURETEX®

V-TPS

FBR

32508 PUR	32508 V-TPS	32508 FBR	8
32510 PUR	32510 V-TPS	32510 FBR	10
32512 PUR	32512 V-TPS	32512 FBR	12
32514 PUR	32514 V-TPS	32514 FBR	14
32516 PUR	32516 V-TPS	32516 FBR	16



Ø 3.75 mm PURETEX®

V-TPS

FBR

37508 PUR	37508 V-TPS	37508 FBR	8
37510 PUR	37510 V-TPS	37510 FBR	10
37512 PUR	37512 V-TPS	37512 FBR	12
37514 PUR	37514 V-TPS	37514 FBR	14
37516 PUR	37516 V-TPS	37516 FBR	16



Ø 4.00 mm PURETEX®

V-TPS

FBR

40008 PUR	40008 V-TPS	40008 FBR	8
40010 PUR	40010 V-TPS	40010 FBR	10
40012 PUR	40012 V-TPS	40012 FBR	12
40014 PUR	40014 V-TPS	40014 FBR	14
40016 PUR	40016 V-TPS	40016 FBR	16



Ø 4.90 mm PURETEX®

V-TPS

FBR

49008 PUR	49008 V-TPS	49008 FBR	8
49010 PUR	49010 V-TPS	49010 FBR	10
49012 PUR	49012 V-TPS	49012 FBR	12
49014 PUR	49014 V-TPS	49014 FBR	14




Ø 6.50 mm

V-TPS


FBR

65010 V-TPS	65010 FBR	10
65012 V-TPS	65012 FBR	12
65014 V-TPS	65014 FBR	14

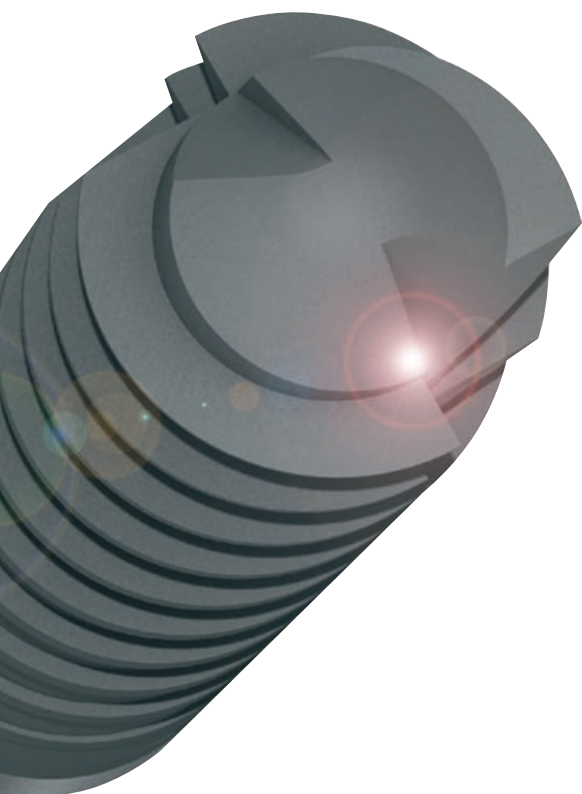

| Simultaneous Drills

Simultaneous drills with fixed stop serve to precisely prepare the implant bed. They are color coded on the shaft center part and indicate the length of the cutting depth. In order to guarantee the optimal osseointegration, they perform a standardized bore shaft adjusted to the length and diameter of the corresponding implant body. The simultaneous drills are slightly undersized, enabling the self-tapping property of the Pitt-Easy implant concept to become effective. The thread profile of the drills favors the conserving of vital bone. Simultaneous drills have an additional length of 1 mm above the graduation-mark as compared to the corresponding implant.

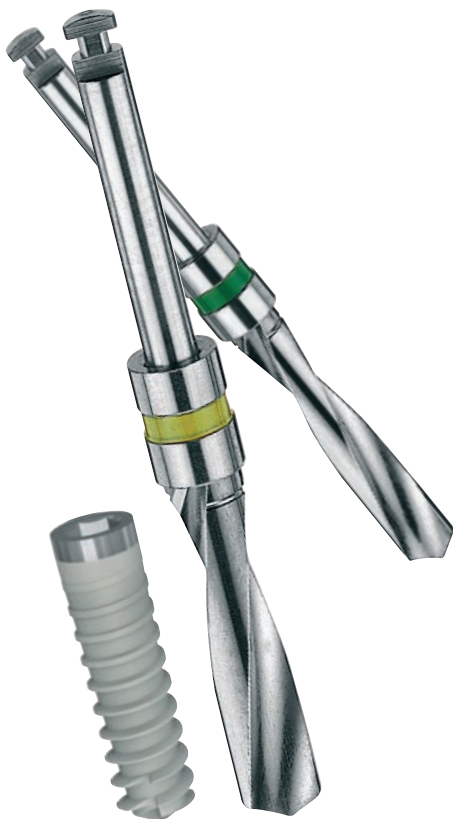

Simultaneous Drills with fixed stop

Item No. Cutting Depth mm

Ø 3.25 mm

92508	8
92510	10
92512	12
92514	14
92516	16




Ø 3.75 mm

97508	8
97510	10
97512	12
97514	14
97516	16



Ø 4.00 mm

94008	8
94010	10
94012	12
94014	14

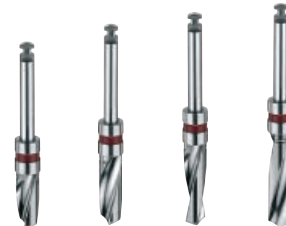

Ø 4.90 mm

94908	8
94910	10
94912	12
94914	14



Ø 6.50 mm

Drills for Ø 6.50 mm see pages 11/12.


| Initial Drills

Initial Drills for all diameters	Item No.	Ø mm	Cutting Depth mm	
Crestotom for the removal of soft tissue, for smoothing and levelling of the alveolar crest	1452			
Drill Extension with latch shaft	1480			
Initial Perforation (IP) Drill for initial perforation of the compact cortical bone	90024	Ø 2.00	16	
Twist Drill, short laser marked 8/10/12/14/16 mm	90422	Ø 2.00	16	
2-Caliber Drill to enlarge the osteotomy from Ø 2.00 mm to Ø 3.00 mm	90413		11	



| Universal Drills

Final Drills	Item No.	Ø mm	Cutting Depth mm	
Final Drill for Ø 3.25 mm one green ring laser marked: 8/10/12/14/16 mm	90426	Ø 3.15 	16	
Final Drill for Ø 3.75 mm one yellow ring laser marked: 8/10/12/14/16 mm	90430	Ø 3.55 	16	
Final Drill for Ø 4.00 mm one red ring laser marked: 8/10/12/14/16 mm	90418	Ø 3.80 	16	
Final Drill for Ø 4.90 mm one blue ring laser marked: 8/10/12/14 mm	90423	Ø 4.60 	14	
Final Drill for Ø 6.50 mm one grey ring laser marked: 8/10/12/14 mm	90446	Ø 6.20 	14	

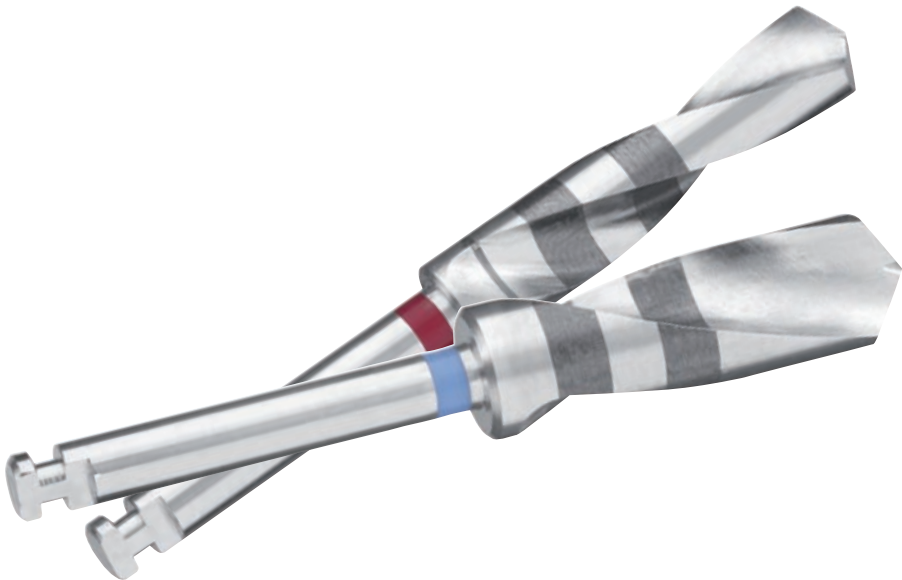

| Universal Drills

Interspace Drills	Item No.	Ø mm	Cutting Depth mm	
Interspace Drill one ring blue/one ring white as intermediate size between 90418 (pg.11) and 90423 (pg.11) laser marked: 8/10/12/14 mm (only for implants-Ø 4.90 mm)	90428	Ø 4.25 	14	
Interspace Drill one ring grey/one ring white as intermediate size between 90423 (pg.11) and 90448 (pg.12) laser marked: 8/10/12/14 mm (only for implants-Ø 6.50 mm)	90447	Ø 5.10 	14	
Interspace Drill one ring grey/two rings white as intermediate size between 90447 (pg.12) and 90446 (pg.11) laser marked: 8/10/12/14 mm (only for implants-Ø 6.50 mm)	90448	Ø 5.60 	14	

12

Cortical Drills	Item No.	Ø mm	Cutting Depth mm	
for compact cortical bone in the mandible				
Cortical Drill for Ø 3.25 mm two green rings	90095	Ø 3.20 	4	
Cortical Drill for Ø 3.75 mm two yellow rings	90096	Ø 3.65 	4	
Cortical Drill for Ø 4.00 mm two red rings	90097	Ø 3.90 	4	
Cortical Drill for Ø 4.90 mm two blue rings	90098	Ø 4.70 	4	
Cortical Drill for Ø 6.50 mm two grey rings	90099	Ø 6.35 	4	

| Surgical Trays

**The new
small tray**



The surgical tray for Ø 3.25 – 4.90 mm

As an alternative to the large surgical tray containing the complete set of simultaneous drills, a small version of the tray is now available for the implant diameters 3.25/3.75/4.00 and 4.90 mm. The small tray contains all components which are required for a quick and rational insertion. An illustrated guide for the assistant staff members is contained in the surgical tray.

Silicone inserts secure the drills and instruments and avoid loss.

Tray Size: 195 mm x 147 mm x 45 mm

The tray contains the following products in insertion sequence to ensure a rational smooth surgery:

- 1 Initial Perforation Drill
- 1 Twist Drill, short
- 1 2-Caliber Drill
- 1 Depth Gauge
- 4 Final Drills and 4 Cortical Drills for Ø 3.25/3.75/4.00 and 4.90 mm
- 1 Interspace Drill for Ø 4.90 mm
- 4 Implant Drivers for Finger Keys or Torque Wrench 2 long and 2 short for Ø 3.25 – 4.90 mm
- 3 Direction Indicators
- 2 Hex Drivers, shaft length 6 and 12 mm
- 1 Drill Extension
- 1 Finger Key Ø 20 mm, for Hexagonal Insert
- 2 O-rings for Prosthetic Drivers
- 1 Driver Sleeve Ø 10 mm
- 1 Pitt-Easy X-Ray Indicator

| Surgical Trays

**For all five
diameters**

Surgical Tray for five diameters

This surgical tray for all five implant diameters offers several improvements as compared to the former version.

The insertion instruments were designed considering the dentists experience and wishes, the trays were equipped accordingly. The tray is delivered with drills and instruments for the diameters 3.25 to 4.90 mm.

Silicone inserts secure the drills and instruments and avoid loss.

An additional drill set for the diameter 6.50 mm is available (see page 15).

The tray contains the following products in insertion sequence to ensure a rational smooth surgery:

- 1 Crestotom
- 1 Depth Gauge
- 3 Initial Drills
- 18 Simultaneous Drills
- 4 Final Drills and 4 Cortical Drills for the diameters 3.25 to 4.90 mm
- 1 Interspace Drill for diameter 4.90 mm
- 4 Finger Keys, manual use
- 2 Implant Drivers for Finger Keys or Torque Wrench
- 6 Parallel Indicators
- 1 Drill Extension
- 1 Hex Driver, driver width 1.70 mm, for manual use
- 4 free insert spaces for drills for diameter 6.50 mm
- 1 Pitt-Easy X-Ray Indicator
- 1 Driver Sleeve

| Surgical Trays

Surgical Trays	Item No.	Ø mm
Pitt-Easy Surgical Tray large surgical tray for all 5 diameters, equipped for Ø 3.25 – 4.90 mm	90212	Ø 3.25  Ø 3.75  Ø 4.00  Ø 4.90 
Surgical Tray, empty	90211	
Set of additional drills for Ø 6.50 mm Final Drill (90446) Interspace Drill (90447) Interspace Drill (90448) Cortical Drill (90099)	90200	Ø 6.50 
Spare filters for large surgical tray (90212), 25 pcs.	90202	
Pitt-Easy Surgical Tray equipped for Ø 3.25 – 4.90 mm	90208	Ø 3.25  Ø 3.75  Ø 4.00  Ø 4.90 
Surgical Tray for Ø 3.25 – 4.90 mm empty	90207	



| Instruments and Accessories

Finger Keys, one-piece

Item No.

Ø mm

Width/Length mm

Finger Key

with rotation control mark and fixed hexagonal insert for manual implant insertion

Finger Key short

for implants Ø 3.25 and 3.75 mm

90093

Ø 3.25

Ø 3.75

2.20

8

for implants Ø 4.00, 4.90 and 6.50 mm

90094

Ø 4.00

Ø 4.90

Ø 6.50

2.50

8



Finger Key long

for implants Ø 3.25 and 3.75 mm

90053

Ø 3.25

Ø 3.75

2.20

16

for implants Ø 4.00, 4.90 and 6.50 mm

90054

Ø 4.00

Ø 4.90

Ø 6.50

2.50

16



Implant Driver

Item No.

Ø mm

Width/Length mm

Implant Driver

for torque wrench (90239, pg.17) – **Please note:**

Strictly follow the insertion instructions.

or Finger Key (90440, pg.17)

Implant Driver short

for implants Ø 3.25 and 3.75 mm

90433

Ø 3.25

Ø 3.75

2.20

8

for implants Ø 4.00, 4.90 and 6.50 mm

90434


Ø 4.00

Ø 4.90

Ø 6.50

2.50

8

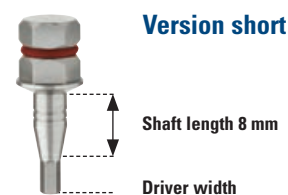

Implant Driver long

for implants Ø 3.25 and 3.75 mm

90431

Ø 3.25

Ø 3.75

2.20

16

for implants Ø 4.00, 4.90 and 6.50 mm

90432


Ø 4.00

Ø 4.90


Ø 6.50


2.50

16

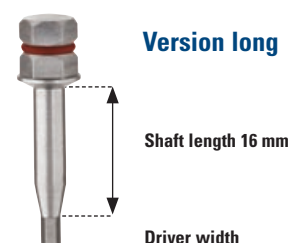

| Instruments and Accessories

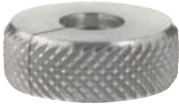	Item No.	Ø mm	
Torque Wrench incl. adjustment screw, with adjustable torque (10, 15, 20, 25, 30, ∞ Ncm)	90239		
	90440		
Implant Driver, Latch-Type for Contra Angle for Ø 3.25 and 3.75 mm 2.2 mm width for Ø 4.00, 4.90 and 6.50 mm 2.5 mm width	90331	Ø 3.25  Ø 3.75 	
	90332	Ø 4.00  Ø 4.90  Ø 6.50 	
Depth Gauge for checking the insertion depth, laser marked: 8 – 16 mm, Ø 2.00 mm	90540		
Direction Indicator center ring Ø 3.50 mm, Ø 2.00 and 3.00 mm	90039		

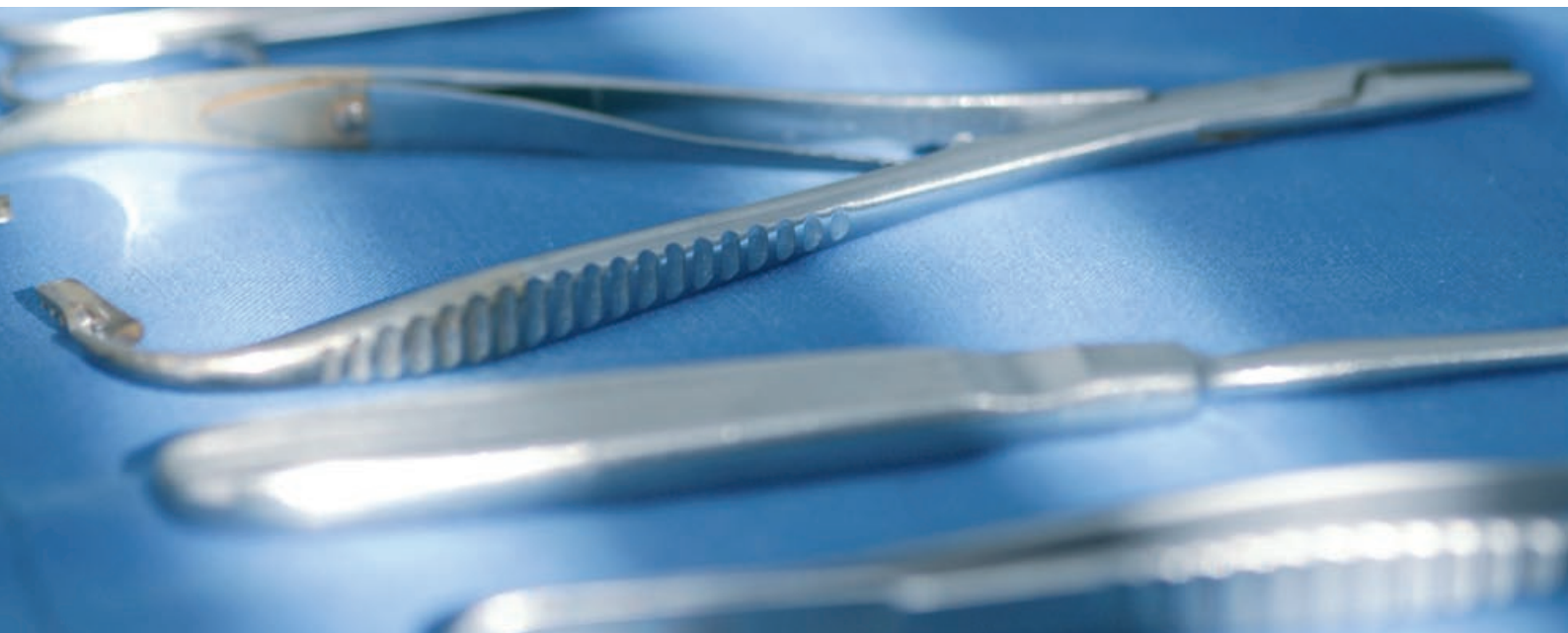

Prosthetics



Open or closed impression

Owing to its flattened side and its horizontal and vertical grooves (version for closed impression) the titanium impression coping is secure against rotation and guarantees optimal repositioning.

In the case of a subgingival implant restoration, the healing abutment can be used to form the tissue. If a two-stage surgical option is chosen, the implant can be uncovered prior to the end of the healing process and the healing abutment can be inserted. After two or more weeks, the soft tissue is healed and an impression can be taken.

For impression taking, the healing abutment is removed and the impression coping is screwed into the implant with the hex driver.

To facilitate repositioning in the impression and to permit control by the laboratory technician after fabrication of the model, the *flat* side of the impression coping should be placed into the implant labially or vestibularly before the impression is made.

For taking impressions only individual impression trays should be used.

The decision as to whether a closed or an open impression method should be used is dependent on the parallelism, number, and angulation of the implants. In the case of a pronounced lack of parallelism (also in relation to the remaining teeth) an open impression is recommended.

The use of healing abutments with conical shape is recommended for prosthetic reconstructions with:

- Locator System
- Paracentric Line
- Ball Head Abutments




Hex Driver with
Driver Sleeve

Impression coping
with implant analog
for closed impression


Impression coping
with implant analog
for open impression




Healing Abutments for gingiva thicknesses from 2 to 6 mm

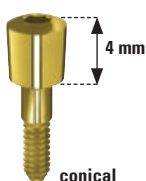


conical

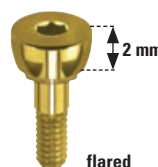

flared

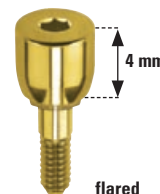

flared

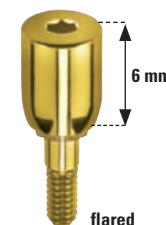

flared

Impressions	Item No.	Ø mm
Impression Coping for open impression including long impression pin	90155	Ø 3.25 
	90156	Ø 3.75 
	90157	Ø 4.00 
	90158	Ø 4.90 
	90159	Ø 6.50 
Impression Coping for closed impression including short impression pin	90055	Ø 3.25 
	90056	Ø 3.75 
	90057	Ø 4.00 
	90058	Ø 4.90 
	90059	Ø 6.50 
Implant Analog	90065	Ø 3.25 
	90066	Ø 3.75 
	90067	Ø 4.00 
	90068	Ø 4.90 
	90069	Ø 6.50 

Hex Driver	Item No.	Length mm
for torque wrench	} use with driver sleeve or torque wrench (pg. 42...)	90625 24
for torque wrench		90623 31
for contra angle (latch)		90323

Healing Abutments	Item No.	Ø mm
Healing Abutment flared gingiva height: 2 mm	32541	Ø 3.25 
	37541	Ø 3.75 
	40041	Ø 4.00 
	49041	Ø 4.90 
	65041	Ø 6.50 
Healing Abutment flared gingiva height: 4 mm	32543	Ø 3.25 
	37543	Ø 3.75 
	40043	Ø 4.00 
	49043	Ø 4.90 
	65043	Ø 6.50 
Healing Abutment flared gingiva height: 6 mm	32546	Ø 3.25 
	37546	Ø 3.75 
	40046	Ø 4.00 
	49046	Ø 4.90 
Healing Abutment conical gingiva height: 4 mm	32544	Ø 3.25 
	37544	Ø 3.75 
	40044	Ø 4.00 
	49044	Ø 4.90 

Organizer

Box with compartments equal to the tooth sequence
inserts (10 pieces each) for spare order

90354

90353



| Abutment Trial Set

Contents of Set	Set-Item No.	Single-Item No.
-----------------	--------------	-----------------

Abutment Trial Set, complete 90347
for checking the optimal position of the abutments.
The set contains 4 trial abutments each
as follows:

- | | |
|--|-------|
| • Trial Abutment straight for Ø 3.25/3.75 mm | 90368 |
| • Trial Abutment straight for Ø 4.00/4.90/6.50 mm | 90369 |
| • Trial Abutment 15° Version A for Ø 3.25/3.75 mm | 90360 |
| • Trial Abutment 15° Version A for Ø 4.00/4.90/6.50 mm | 90361 |
| • Trial Abutment 15° Version B for Ø 3.25/3.75 mm | 90362 |
| • Trial Abutment 15° Version B for Ø 4.00/4.90/6.50 mm | 90363 |
| • Trial Abutment 25° Version A for Ø 3.25/3.75 mm | 90364 |
| • Trial Abutment 25° Version A for Ø 4.00/4.90/6.50 mm | 90365 |
| • Trial Abutment 25° Version B for Ø 3.25/3.75 mm | 90366 |
| • Trial Abutment 25° Version B for Ø 4.00/4.90/6.50 mm | 90367 |

Trial abutments are also available separately as
set of four (see single-item no.)

Detailed description with cleaning instructions
is included.



The Abutment Trial Set serves to predetermine the exact prosthetic abutments and positioning in practice and laboratory.



| A.G.T.® Line

Esthetic Temporary Abutment

The A.G.T.® abutment is used for post-operative gingival moulding, and at the same time for initiating esthetic temporary treatment. Its design is preconstructed with an anatomical gingival shoulder. The abutment is made of a tooth-coloured acrylic material and, like the final titanium or ceramic abutment which can be selected later, can be fastened precisely with a retaining screw to prevent rotation. It is used for the fabrication of a temporary crown and should only remain in the mouth until the gingiva has healed.

Temporary cement can be used for cementation.

After healing and formation of the papillae, the temporary crown is easily removed by loosening the retaining screw along with the A.G.T.® abutment.



Impression – see page 18

A.G.T.® Abutment

Item No.

Ø mm

A.G.T.® Abutment
acrylic (P.O.M.)

32588

Ø 3.25

37588

Ø 3.75

40088

Ø 4.00

49088

Ø 4.90

Hex Driver

Item No.

driver width 1.70 mm with friction grip,
for retaining screw,
use with driver sleeve
or torque wrench (pg. 42...)

90623

31

21

Retaining screw
can be shortened
up to the marking



| V.D.L. Esthetic Line


The titanium anti-rotation abutments of the V.D.L. Esthetic Line show several advantages. The most important feature is the circular, groovelike shoulder which is lower in the buccal area for esthetics and slightly ascends with an elegant swing towards the lingual. This enables the technician to prepare an exact, clearly defined crown margin. The subgingival design of the V.D.L. Esthetic abutments allows individual adjustment to the gingiva contour as well as an adjustment to the anatomical direction of the tooth.

The internal design of the abutment was advanced so that the outer shape of the abutment remains slim, which is of considerable advantage for the design of narrow incisors. We recommend the use of special titanium drills for the trimming of the V.D.L. Esthetic abutments.


The retaining screw is secured by an additional internal thread system to avoid losing the screw or mistaking it for another. For cleaning purposes, the retaining screw can be unscrewed from the abutment, if necessary.

Please note: At the time of the implant insertion the abutment position of Version A or B should be taken into consideration (flat side/version A or corner/version B showing towards labial/buccal). The final trial is made with the original abutment or with a trial abutment (see page 20).

The design of the V.D.L. Esthetic abutment enables the technician to meet and exceed the high demands of function and esthetics of today's implant prosthetics.

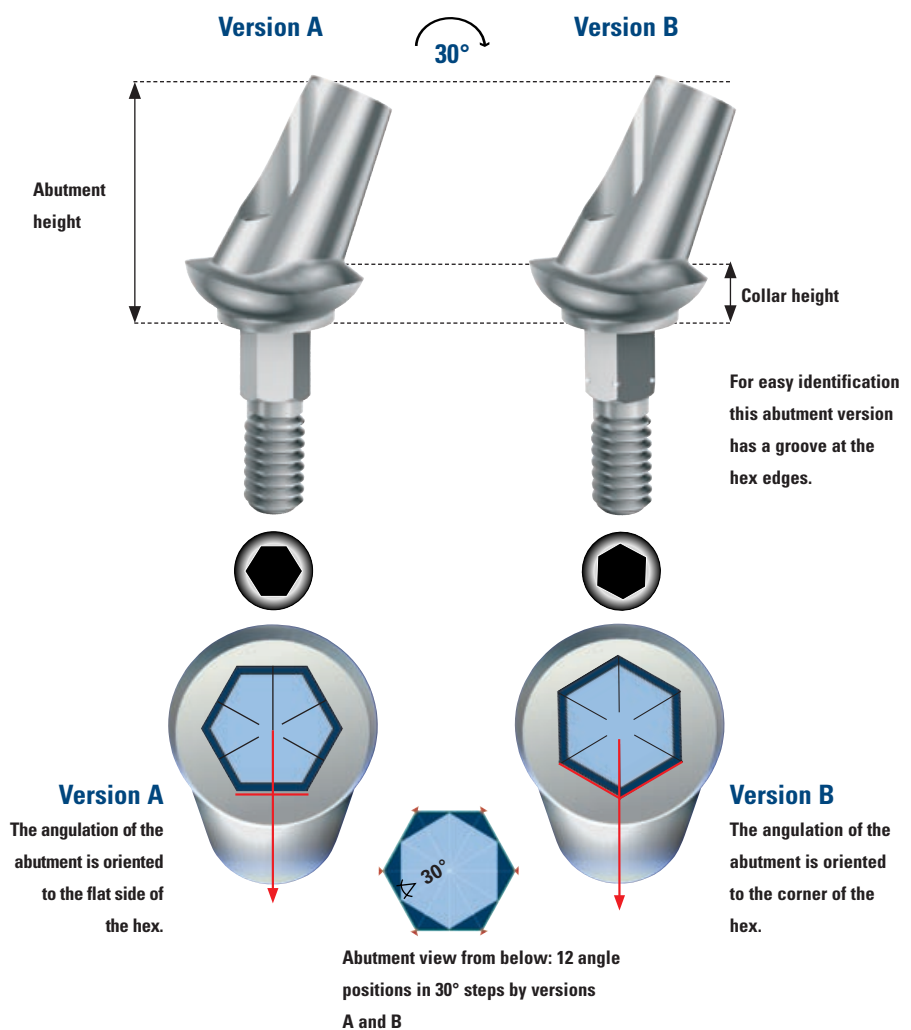

Impression – see page 18

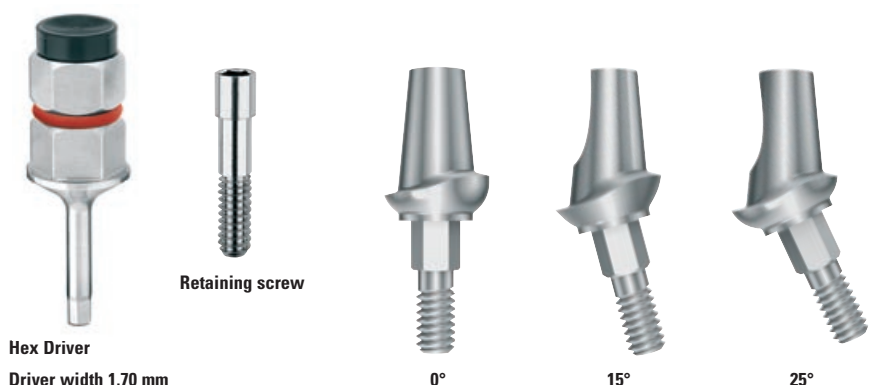

| V.D.L. Esthetic Line

V.D.L. Esthetic Abutments	Item No.	Item No.	Ø mm
V.D.L. Esthetic Abutment, titanium, 0° straight collar height 1.20 mm/abutment height 7.50 mm	 = A	 = B	
	32656		Ø 3.25 
	37656		Ø 3.75 
	40156		Ø 4.00 
	49156		Ø 4.90 
	65156		Ø 6.50 
V.D.L. Esthetic Abutment, titanium, 0° straight collar height 2.00 mm/abutment height 8.30 mm	32659		Ø 3.25 
	37659		Ø 3.75 
	40159		Ø 4.00 
	49159		Ø 4.90 
V.D.L. Esthetic Abutment, titanium, 15° angled collar height 1.20 mm/abutment height 8.20 mm	32654	32675	Ø 3.25 
	37654	37675	Ø 3.75 
	40154	40175	Ø 4.00 
	49154	49175	Ø 4.90 
	65154		Ø 6.50 
V.D.L. Esthetic Abutment, titanium, 25° angled collar height 1.20 mm/abutment height 8.20 mm	32655	32670	Ø 3.25 
	37655	37670	Ø 3.75 
	40155	40170	Ø 4.00 
	49155	49170	Ø 4.90 
Hex Driver	Item No.	Length mm	
driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...) for contra angle (latch)	90623	31	
	90323		

| Castable Abutment System

The C.A.S. Novobase abutment and the C.A.S. Esthetic abutment can be used for many different indications for dentist-removable or cemented superstructures. They allow an individual design of the implant reconstruction which corresponds with the individual patient situation. Even extremely unfavorable angulations or axial divergencies (that are beyond the possibility of confectioned abutments) can be compensated by using this system. Also for telescopic or conus crowns the C.A.S. abutments can be used for perfect results.

The abutments can be shortened, individually modelled and also angulated.

For smoothing the internal screw channel should be prepared with a reamer (Item No. 90247).

C.A.S. Novobase® for a perfect accuracy

The C.A.S. Novobase abutment is a cast-on version. The base consists of machined gold-platinum, so that all precious metal alloys can be used for a cast-on method. The CNC-drilled connection part to the implant guarantees an absolutely precise fit. The upper part is made of a castable acrylic. For direct porcelain application, the use of a ceramic alloy is essential. The abutments can be shortened, individually modelled and also angulated. For an exact fit of the retaining screw, the internal screw channel should be prepared with a reamer (Item No. 90247, pg. 25).

Note: If the complete height of the shaft is used, the long hex driver should be used (Item No. 90623). For a lower height – below the circular groove marking – a short hex driver is sufficient (Item No. 90625)

Note for processing

Interval for fusing: 1,400 – 1,460° C
Preheating temperature: up to 850° C
Casting temperature: up to 1,200° C

C.A.S. Novobase® Abutments	Item No.	gram*	Ø mm
C.A.S. Novobase® Abutment acrylic/precious metal (Novostil)	32527	0.43	Ø 3.25
	37527	0.47	Ø 3.75
	40027	0.55	Ø 4.00
	49027	0.63	Ø 4.90

*precious metal component in gram

Hex Driver	Item No.	Length mm
driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...)	90623	31

Retaining screw



Acrylic, castable



Precious metal

The investment material should be removed by using blast brightening pearls at low blasting pressure of ≤ 1 bar.



Impression – see page 18

| Castable Abutment System

C.A.S. Esthetic

The C.A.S. Esthetic abutment of castable acrylic can be individually trimmed, adjusted and waxed up in any dimension on the model prior to casting. Also the angulation can be changed by wax-up.

The C.A.S. Esthetic abutment is castable with all precious metal alloys.

The design of this acrylic abutment with anatomical shoulder follows the V.D.L. Esthetic Line of titanium (pg. 22/23). It can easily be adjusted to the present marginal gingiva contour after gingiva healing.

The shape of this abutment eliminates the necessity to perform a wax-up or other modifications in many cases. This facilitates the handling for the technician and saves time. The abutment can be used for all casting alloys (titanium only by pressure casting procedure).

We recommend the use of precious metal alloy!

The preparation of casting channels should always be started laterally in the direction of the hex. The hex guarantees a perfect anti-rotation safety and is able to shape itself after casting and thus ensure stability. Please pay attention to exact modelling and trimming of the hex in order to obtain a perfect fit to the implant. For the precise fit of the retaining screw, the internal screw channel should be prepared with the reamer. For fixation and unscrewing of the retaining screw, please use the hex driver.

C.A.S. Esthetic Abutments

Item No.

Ø mm

C.A.S. Esthetic Abutment

acrylic, castable

32537

Ø 3.25

37537

Ø 3.75

40037

Ø 4.00

49037

Ø 4.90

Reamer

Item No.

for the Castable Abutment System and the Paracentric Bar Element Acrylic for finishing of the screw channel after casting

90247



Hex Driver

Item No.

Length mm

driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...)

90623

31

Retaining screw can be shortened up to the marking



CAD/CAM-Base

New

CAD/CAM-Base

The CAD/CAM-Base abutment serves to fabricate individual zirconium abutments for the Pitt-Easy implant system.

This prosthetic construction puts the user in a position to satisfy the highest esthetic demands placed on implant-supported dentures. The abutment is constructed in two parts so that in production the fabricated zirconium part is bonded with the titanium base.

To make an individual wax-up easier, a plastic chimney is placed on the titanium base as a modelling base at the outset.

In order to represent the subsequent screw channel, a scan screw is optionally available, which can be screwed with the titanium base for scanning the "inner contours".

Important:

The instructions must be read completely prior to processing the product!

CAD/CAM-Base	Item No.	Ø mm
--------------	----------	------

CAD/CAM-Base	32569	Ø 3.25
	37569	Ø 3.75
	40069	Ø 4.00
	49069	Ø 4.90

supplied incl. laboratory screw and retaining screw

Screws	Item No.	Length mm
--------	----------	-----------

Retaining screw	90276	10 mm
Laboratory screw	90278	10 mm
Scan screw	90282	17 mm

Instruments	Item No.
-------------	----------


Hex Driver

for torque wrench, shaft length 6 mm
for torque wrench, shaft length 12 mm
for contra-angle

90625
90623
90323

Torque Wrench

90239

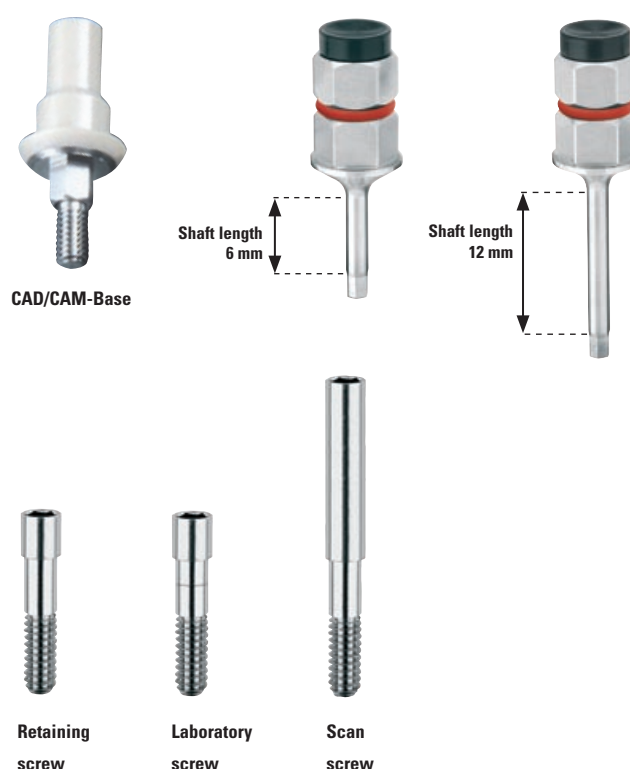



| Individual Abutments

With the Individual Abutments we have created abutments for which the primary parts for the telescopic or conus crown technique can be drilled from the full material. All implant analogs can be provided with the abutments corresponding to the diameters.

The Individual Abutments are ideally suitable also for the preparation of individual angulated abutments, for instance if a abutment is required which should have another angulation than 0°, 15° or 25°. Also, in case a special gingiva contour requires a circular shoulder which cannot be prepared from a confectioned abutment, this solid Individual Abutment can be used. The crown or bridge construction can be cemented or individually screwed to this Individual Abutment.

By means of direct veneering with titanium ceramic material individual ceramic abutments or directly screw-retained crowns can be fabricated.

For trimming of the Individual Abutment on the model or with the H.U.G.O. instrument (page 29), the laboratory screws are provided. These have a circular marking at the screw shaft and should be used during the lab procedure. The retaining screws are used for final fixation of the prosthetic construction in the mouth. We recommend the use of titanium drills.

When planning a telescoping construction on implants or implants in combination with natural teeth, please pay attention to the construction's static safety, i.e. distal support.

Laboratory screw



Retaining screw



Individual Abutment, 11 mm



Individual Abutment, 7.5 mm



Retaining screw inserted

Individual Abutment, inserted in H.U.G.O instrument



Individual Abutment, trimmed



Impression – see page 18

| Individual Abutments

Individual Abutments	Item No.	Ø mm
Individual Abutment short length 7.5 mm, titanium, incl. retaining screw and lab screw	32549	Ø 3.25 
	37549	Ø 3.75 
	40049	Ø 4.00 
	49049	Ø 4.90 
	65049	Ø 6.50 
Individual Abutment long length 11 mm, titanium, incl. retaining screw and lab screw	32564	Ø 3,25 
	37564	Ø 3,75 
	40064	Ø 4,00 
	49064	Ø 4,90 

Screws	Item No.	Length mm
Retaining screw	90276	10
Laboratory screw	90278	10

Hex Driver	Item No.	Length mm
driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...)	90623	31
for contra angle (latch)	90323	

Set of instruments	Item No.
H.U.G.O. set of holding instruments for abutments holding, grinding, smoothing, optimizing	90250



| Paracentric Line

Paracentric Line is especially designed for bar structures with diverging abutments.

This unique system provides the opportunity to use screw-retained superstructures with easy compensation of divergencies. The Paracentric Line consists of three parts: permucosal extension, bar connection element, and retaining screw. The titanium permucosal extension with hexagonal base fits precisely into the implant body's hexagonal socket. The

different heights of the permucosal extensions compensate the varying thickness of the mucosa. The bar element is available either in acrylic, precious metal, castable alloy, or titanium. Since the hex provides primary stability to the permucosal extension and the bar element is being stabilised through the cone, the forces applied to the bar construction are diverted from the retaining screw and prevent loosening and fracture.



Hex Driver

Retaining screw
Titanium

Bar Element



Acrylic



Titanium



Precious Metal

Advantages for the Dentist

- easy removal of bar for proper cleaning at any time
- permucosal extension and bar retained by one retaining screw
- excellent superstructure evaluation at try-in due to horizontal plane above gingival level (gap detection)

Advantages for the Laboratory Technician

- the 20° cone angle permits a corresponding correction of divergences
- due to the different materials of the bar segment the technician can cast the construction (acrylic precast parts); cast on, solder, or glue semiprecast parts and laser-weld or glue semiprecast titanium components

Note for processing

Interval for fusing: 1,400 – 1,460° C
Preheating temperature: up to 950° C
Casting temperature: up to 1,200° C



Impression – see page 18



Bar Element

Permucosal
Extension
TitaniumImplant
analog

≤ 40°


Implant divergencies up to 40° can be compensated by the conical design of the connecting piece of the permucosal extension to the bar segment.

| Paracentric Line


Paracentric Abutments	Item No.	Ø mm
Paracentric Bar Element		
5 mm height		
acrylic	90277	
precious metal (Novostil)	90279	
titanium (not recommended for titanium casting)	90281	
Paracentric Permucosal Extension		
titanium	32577	Ø 3.25 
	37577	Ø 3.75 
with retaining screw, length 14 mm	40077	Ø 4.00 
2 mm height	49077	Ø 4.90 
Paracentric Permucosal Extension		
titanium	32578	Ø 3.25 
	37578	Ø 3.75 
with retaining screw, length 16 mm	40078	Ø 4.00 
4 mm height	49078	Ø 4.90 
Paracentric Permucosal Extension		
titanium	32567	Ø 3.25 
	37567	Ø 3.75 
with retaining screw, length 18 mm	40067	Ø 4.00 
6 mm height	49067	Ø 4.90 
Hex Driver	Item No.	Length mm
driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...)	90623	31
Reamer	Item No.	
for polishing the cast after veneering (for Paracentric Bar Element acrylic)	90247	




| Ball Head Abutment System

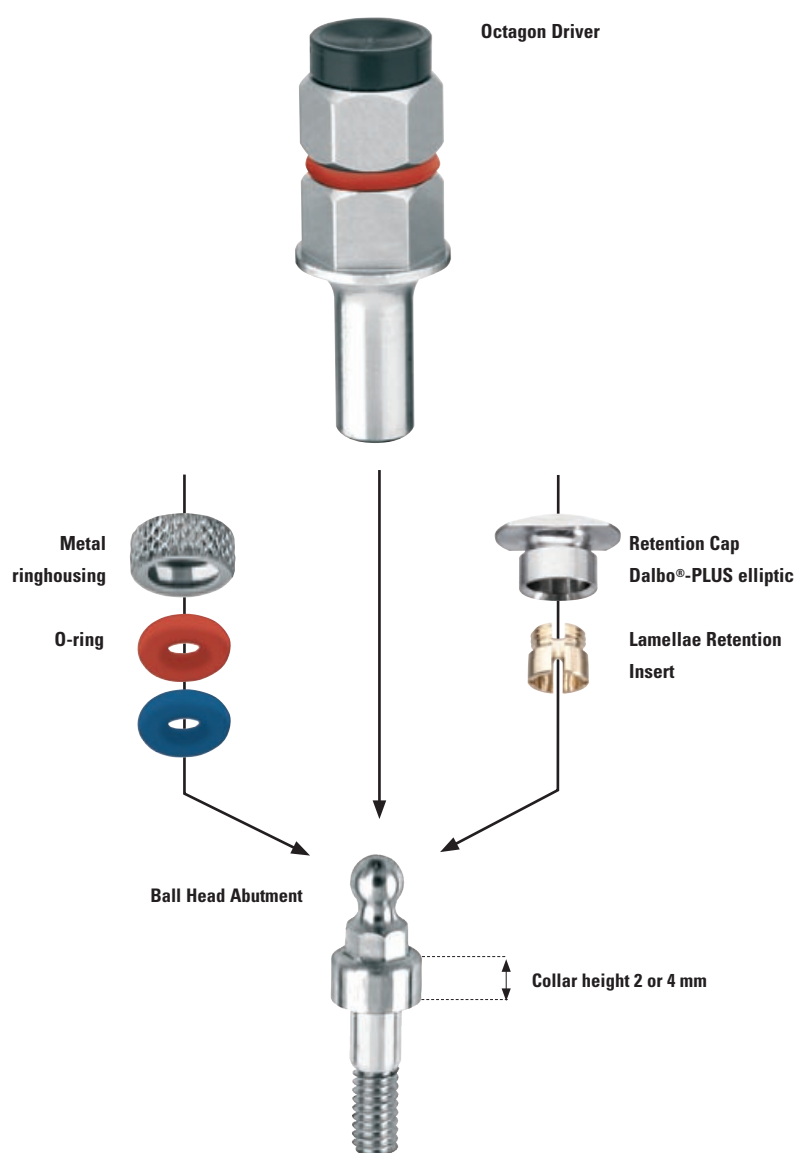

The Ball Head Abutment (titanium) supports a resilient anchoring element, secures the denture from different forces (such as pulling, pressure or lateral load), and guarantees a relief of the implants when rotation forces in different directions are stressing the prosthetic construction.

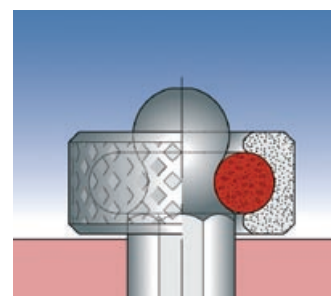

**Metal
Ringhousing
with O-ring**

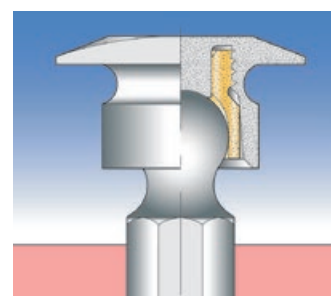

**Retention Cap
Dalbo®-PLUS
elliptic
Shear strength
can be adjusted
from 200 to
1,200 grams**

Indication

- stabilization of a removable prosthesis (preferably supported by four implants in the interforaminal region of the mandible)

Ball Head Abutments are not indicated for connected constructions.

From now on, you have the possibility to use Ball Head Abutments either with the O-ring and metal-ringhousing or with the retention cap Dalbo-PLUS elliptic. For smaller interdental spaces this abutment variation offers an economic and esthetic solution (no gingiva metal transparency).

| Ball Head Abutment System

Ball Head Abutments	Collar Height	Item No.	Ø mm	
Ball Head Abutment titanium	2 mm	32551	Ø 3.25	
	2 mm	37551	Ø 3.75	
	2 mm	40051	Ø 4.00	
	2 mm	49051	Ø 4.90	
	4 mm	32553	Ø 3.25	
	4 mm	37553	Ø 3.75	
	4 mm	40053	Ø 4.00	
	4 mm	49053	Ø 4.90	

Single items	Item No.	Item No.
Metal Ringhousing and O-rings complete	90301	Implant Analog for Ball Head Abutments 2251
Retention Cap Dalbo®-PLUS elliptic* complete, incl. Lamellae Retention Insert	90289	Octagon Key for Ball Head Abutments use with driver sleeve or torque wrench (pg. 42...) 90650
Metal Ringhousing individual	3582	Activator* for retention cap Dalbo®-PLUS elliptic 90288
O-ring individual red 3587 blue 3589		
Lamellae Retention Insert* individual	90287	 Impression – see page 18

| Locator Abutment System

Locator Abutments		Collar	Item No.	Ø mm	
The Locator Abutment is an attachment with a self-aligning feature. This feature will aid patients in seating their denture and eliminate additional wear from improper seating.	1 mm	8311	3.25 mm		
	3 mm	8313	3.25 mm		
	1 mm	8321	3.75 mm		
	3 mm	8323	3.75 mm		
1 mm	8331	4.00 mm			
3 mm	8333	4.00 mm			
With the Locator's minimal vertical height, it is ideal where interocclusal space is limited. This abutment can also accommodate up to 40° between implants.					

34

Locator Attachments	Item No.	Locator Attachments	Item No.	
Locator Male Processing Package 2 Sets, each Set include: 1 Processing Replacement Male 8515, 1 Replacement Male 8524/8527/8529, 1 Spacer white	8519	Locator Extended Range Replacement Male green, 4 Pack	8547 	
Locator Replacement Male clear, 4 Pack	8524 	Locator Extra Light Extended Range Male red, 4 Pack	8548 	
Locator Light Retention Replacement Male pink, 4 Pack	8527 	Locator Black Processing Replacement Male black, 4 Pack	8515 	
Locator Extra Light Retention Replacement Male , blue, 4 Pack	8529 			
Locator Parallel Post 4 Pack	8517	 Locator Parallel Post	 Locator Impression Coping	 Female Analog Ø 4 mm
Locator Impression Coping 4 Pack	8505	 ZEST ANCHORS INC. 800-262-2310 25° 20° 15° 10° 0° 10° 15° 20° 25° Angle Measurement Guide		
Locator Female Analog 4 mm Diameter. 4 Pack	8530			
Angle Measurement Guide	9530			
Locator Core Tool	8393			
Locator Key	8317			
Manufacturer: Zest Anchors Inc.				
Torque Wrench	90239			

| Standard Abutments with hex

The standard abutment (titanium) with hex is indicated for cemented crown and bridge restorations.


The standard abutment is retained directly to the implant by the retaining screw.

Impression taking

The impression is taken by using the impression coping.



Impression – see page 18

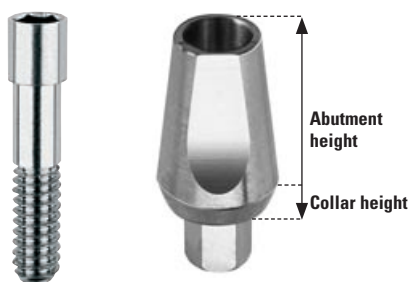

Standard Abutment 0°



Standard Abutment 15°

Standard Abutment 25°

Standard Abutments	Collar/Abutment height	Item No.	Ø mm
Standard Abutment, titanium 0° straight	collar height 1 mm abutment height 7 mm	32556	Ø 3.25
		37556	Ø 3.75
		40056	Ø 4.00
		49056	Ø 4.90
		65056	Ø 6.50
	collar height 3 mm abutment height 9 mm	32559	Ø 3.25
		37559	Ø 3.75
		40059	Ø 4.00
		49059	Ø 4.90
Standard Abutment, titanium 15° angulation Version A (pg. 22)	collar height 1 mm abutment height 9 mm	32554	Ø 3.25
		37554	Ø 3.75
		40054	Ø 4.00
		49054	Ø 4.90
	collar height 3 mm abutment height 11 mm	32557	Ø 3.25
		37557	Ø 3.75
		40057	Ø 4.00
		49057	Ø 4.90
Standard Abutment, titanium 25° angulation Version A (pg. 22)	collar height 1 mm abutment height 9 mm	32555	Ø 3.25
		37555	Ø 3.75
		40055	Ø 4.00
		49055	Ø 4.90

Hex Driver	Item No.	Length mm
driver width 1.70 mm with friction grip, for retaining screw, use with driver sleeve or torque wrench (pg. 42...) for contra angle (latch)	90623	31
	90323	

| Prosthetic Accessories

Bar and Clip System

The Bar and Clip System provides easy external splinting of implants and is particularly suitable for screw-retained bar frameworks and clip-attached prostheses.

The system consists of

- Resilient/Rigid castable Bars
- Clips, made of permanently flexible, high-precision acrylic
- Metal Housings for Clips (chrome-nickel-alloy)
- Press-In Pin

With all accessory parts made of castable acrylic, the laboratory technician is free to choose among all available alloys, ideally titanium, of course. This structure is castable in one piece thus keeping the alloy in the oral cavity at a minimum. The clip housing, made of metal, assures an optimal positioning of the clip.

Metal Housing
for Clips



Clips



Round Bar, 50 mm



Bar, drop profile, 50 mm



Bar, drop profile, Degulor M, 25 mm



Press-In Pin




Profile Bar, rigid, 50 mm



Root-Buffer Attachment

Universally castable ball anchor serves as retention element for all removable implant superstructures, for incorporation into the prosthesis or bar system and can be placed in the center or at end position. Can be included in the wax-up for casting.

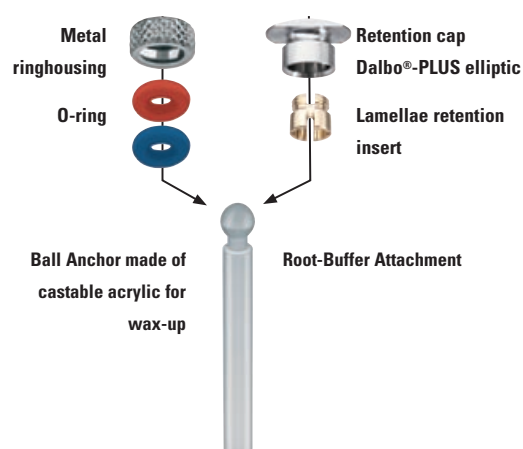



| Prosthetic Accessories

Bar and Clip System	Item No.
Single items, delivery of 6 resp. 4 each	
• 6 Clips, yellow, Ø 1.80 mm (stronger friction)	3577
• 6 Clips, orange, Ø 1.90 mm (weaker friction)	3572
• 4 Metal Housings for Plastic Clips	3566
• 6 Profile Bars, green, L 50 mm	3573
• 6 Bars, drop profile, yellow, L 50 mm	3574
• 6 Round Bars, blue, Ø 1.80 mm, L 50 mm	3575
• 1 Press-In Pin, black	3555
• 1 Bar, drop profile, gold alloy Degulor M	3563
	Item No.
Root-Buffer Attachment	
Ball anchor, castable acrylic, (Ø 2.25 mm)	90298
Metal Ringhousing	3582
O-ring	
red	3587
blue	3589
Retention Cap Dalbo®-PLUS elliptic complete* (for Ø 2.25 mm), producer: Cendres & Métaux SA	90289
Metal Ringhousing and O-rings complete* (for Ø 2.25 mm)	90301

*single items see pg. 33



| Table Prosthetic Lines and Abutments

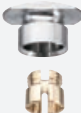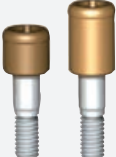I Titanium			
Abutments for impression taking	Abutments without hex for complete denture, not rotation-safe	Abutments with anti-rotational hex	
 Impression Coping closed	 Ball Head Abutment 2/4 mm neck height	 Standard Abutment straight 0°, 15°, 25° anti-rotational 7/9 mm height	
 Impression Coping open	Anchoring elements for Ball Head Abutments  Metal Ringhousing with O-rings	 V.D.L. Esthetic Line straight 0°, 15°, 25° with retention of retaining screw	
 Healing Abutment conical 4 mm height	 Retention cap Dalbo®-PLUS elliptic with lamellae retention insert	 Individual Abutment to be milled individually 7.5 and 11 mm height	
 Healing Abutment flared 2/4/6 mm height		 Paracentric Line Bar Element 5 mm height	
 Implant Analog	 Locator System	 Permucosal Extension 2/4/6 mm height Titanium	

I Precious metal	I Acrylic	I Zircon
Abutments with anti-rotational hex	Abutments with anti-rotational hex	Abutments with anti-rotational hex
 C.A.S. Novobase® Castable Abutment System castable acrylic/precious metal	 C.A.S. Esthetic Castable Abutment System castable	New  CAD/CAM-Base
	 A.G.T.® Line anatomic healing abutment for temporary care	
 Paracentric Line Bar Element 5 mm height	 Paracentric Line Bar Element, castable 5 mm height	
 Permucosal Extension 2/4/6 mm height Titanium	 Permucosal Extension 2/4/6 mm height Titanium	

| Prosthetic Kit



Prosthetic Kit

The basic tray includes a short hex driver (90625) and a driver sleeve (90439). All other components can be ordered separately, as required.

Kit with basic items

90358

Hex Driver short
Driver Sleeve

90625
90439


Additional components (to be ordered individually)

Item No.

Hex Driver long	90623
Hex Driver for contra angle (latch)	90323
Octagon Key	90650
Activator Dalbo Plus	90288
Pitt Locator Key	08317
Hand Grip for technician	90442
Locator Core Tool	08393
Torque Wrench	90239
with adjustment screw 25 Ncm	
H.U.G.O., 2 instruments	90250



Table Prosthetic Drivers and Screws

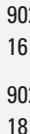I Hex Drivers		I Octagon Key		I Locator Key	
Driver width 1.70 mm  Shaft length 6 mm Shaft length 12 mm					
 Implant Cover Screw	 Healing Abutment	 Impression Coping open	 Impression Coping closed	 Ball Head Abutment	 Locator System
I Screws					
All abutments are delivered incl. screw. For repeat orders please find the corresponding order Nos. in this table.					
 90251 11 mm	 90253 13 mm	 90264 16 mm*	 Scan screw 90282 17 mm	 Retaining screw 90276 10 mm	 Laboratory screw 90278 10 mm
 Standard Abutment with Hex	 A.G.T. Line	 C.A.S. Esthetic	 CAD/CAM-Base	 V.D.L. Esthetic Line	 C.A.S. Novobase
			 Individual Abutment		 Paracentric Line
					 90270 14 mm*
					 90271 16 mm*
					 90272 18 mm*

*can be trimmed

| Prosthetic Drivers, Instruments and Accessories


Item No.

Two-piece driver for safe and easy application

The driver can be used with the torque wrench with adjustable torque or in connection with the driver sleeve if the torque does not have to be adjusted. For final fixation of the abutment in the implant, the torque wrench should be used.

Hex Driver

driver width 1.70 mm for: Implant cover screws, impression abutments, healing abutments and all retaining and laboratory screws

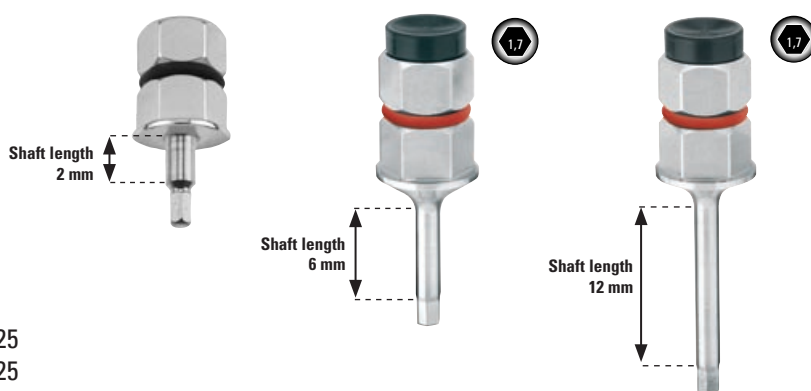

for torque wrench, shaft length 2 mm	90425
for torque wrench, shaft length 6 mm	90625
for torque wrench, shaft length 12 mm	90623
for contra angle	90323

Octagon Key

for Ball Head Abutments
for torque wrench
for contra angle

90650
90324



Locator Key

8317



Driver Sleeve

for manual use of the following prosthetic drivers: 90625, 90623, 90650, 8317.

90439




The driver sleeve can be easily fixed onto the driver.

	Item No.	
Hand Grip for technician In order to provide optimal handling for the technician in the laboratory, a special hand grip was developed. The rotating disk head ensures a relaxed and quick processing. The hand grip fits all drivers which are suitable for torque wrench use.	90442	
Torque Wrench for fixation of the retaining screws (and one-piece-abutments) with adjustable torque (10, 15, 20, 25, 30, ∞ Ncm). For the following prosthetic drivers: 90625, 90623, 90650, 8317. The torque wrench can easily be fixed onto the driver.	90239	
H.U.G.O. set of holding instruments for abutments holding, grinding, smoothing, optimizing	90250	

Instruments and Accessories	Item No.	
Centering Punch Ø 3.50 mm for manual re-entry of the gingiva	1020	
Tissue Punch for reduced contra-angle for exposure of the osseointegrated implants Ø 3.50 mm outer diameter Ø 4.20 mm outer diameter	90085 90086	 
Titanium Tweezers	2852	

PITT-EASY®-Implant

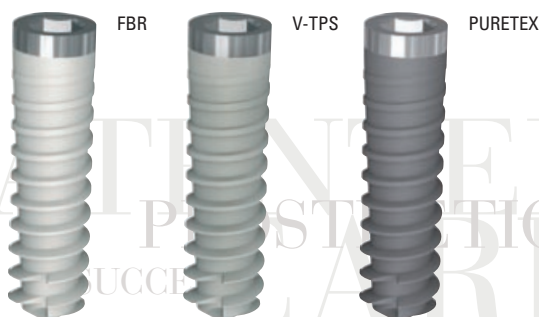

Innovative products

-  PITT-EASY®
-  ENDOPORE®
-  BICORTICAL®
-  BIORESORB®
-  CYTOPLAST®

Sybron Implant Solutions

The dental implant platform of Sybron Dental Specialties (SDS) brings together three (INNOVA – ORALTRONICS – ATTACHMENTS) organizations that have been leaders in implant manufacturing, design and development since the inception of implant dentistry.

Delivering two of the most respected, leading-edge dental implants on the market – ENDOPORE® and PITT-EASY® – we are giving you access to a world of implant solutions.

Sybron Implant Solutions GmbH

Julius-Bamberger-Straße 8a · 28279 Bremen · Germany

Phone +49 421 43939-0 · Fax +49 421 443936

www.sybronimplants.de · info@sybronimplants.de

CE 0482